



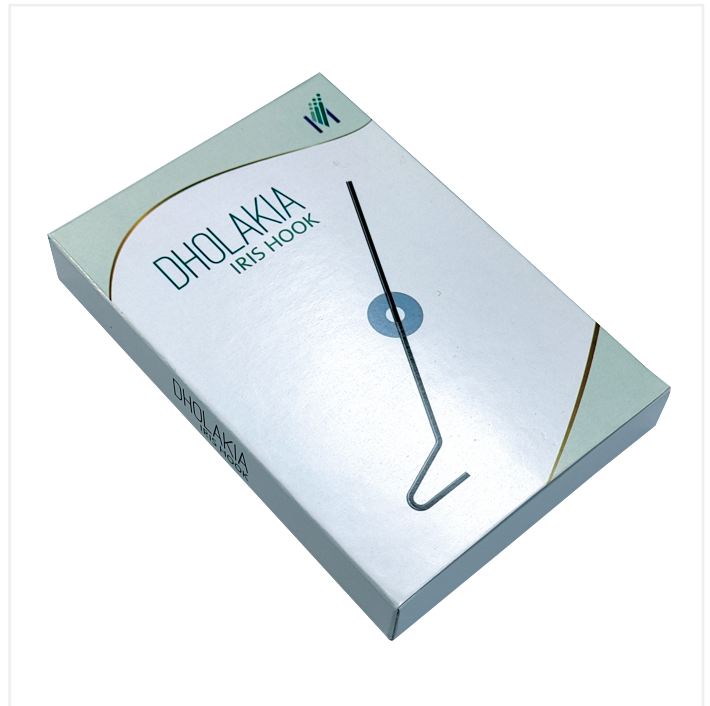

Dholakia Iris Hook

Dholakia Iris Hook is a set of 5 disposable iris retractors made of black nylon fibre with an adjustable silicone lock. The hook has a 'U' shaped bend on one end, which helps in engaging the iris and pulling it backwards to enlarge the pupil. The hook is designed in such a way that it is angled to be in the plane of the anterior capsule. The silicone lock helps in retaining the hook in retracted position during the surgery. After the surgery, the lock can be adjusted to release the iris so that it comes back to its original position.

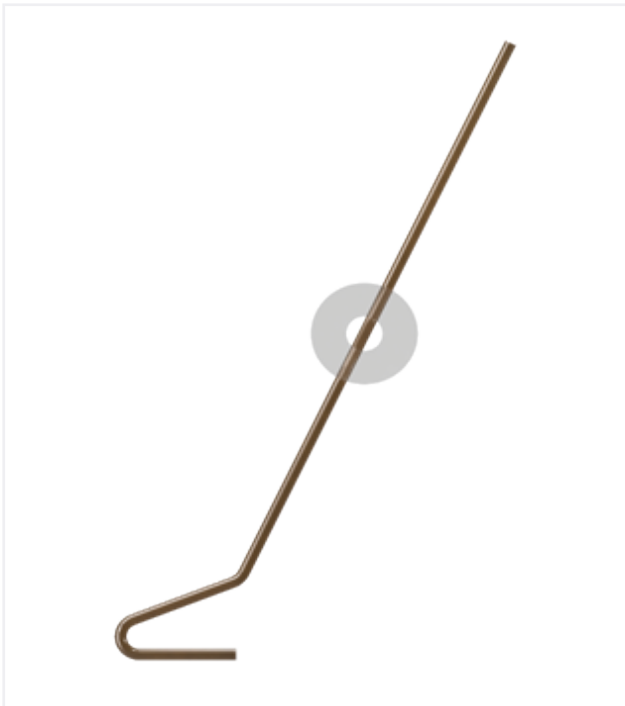

Product Description

Product Code	Description	Qty. per box
MIPL/D16	Dholakia Iris Hook (Nylon & Silicone)	5 Pcs. (Sterile)

Features

- "U" shaped bend to engage the iris.
- Hook lies in the plane of the Anterior Capsule.
- The stem is at acute angle from the plane making it easy to handle.
- Adjustable silicone lock helps in retaining the hooks in a retracted position during the surgery.
- Supplied sterile in a set of 5 hooks in a plastic container inside a medical grade peel open pouch.

Advantages

- Unique angled design diminishes possibilities of the anterior chamber becoming shallow upon iris retraction.
- Reduces the chances of Sphincter damage.
- Makes IOL insertion easy.
- Reduces chances of Endothelium damage.
- Easy to insert and remove.

Indications

Iris Hooks provides temporary retraction (pupil dilation) of the iris where dilation cannot be achieved pharmacologically.

Instructions for Use

The instructions for use and surgical technique are supplied in the form of a leaflet with instructions given in diagram form where appropriate. The international symbols used on the packaging are also explained.

Sterilization

Products are sterilized through Ethylene Oxide (EO) according to a validated cycle.



MADHU INSTRUMENTS PVT LTD

A-260, Okhla Industrial Area, Phase - 1, New Delhi - 110020, India

T: (+91-11) 4270 1030, 31 | E: sales@madhuinstruments.com | www.madhuinstruments.com

Customer Care No. +91 9810866176 (for international queries), +91 7428186111 (for domestic queries)