

Micro Incision Blade with Safety

A specially designed handle with a retractable safety cap that protects the blade and minimises the risk of accidental ophthalmic blade injuries that could occur during all phases of a surgical procedure. A simple slide mechanism pulls the safety cap back making the blade visible and easy to use. The safety cap locks into place during the surgery. Once the surgery is completed the safety cap is pushed forward to bring the blade back under the safety cap.




Blade in Safety Cap



Blade ready to use


Lance Tip Blade These blades are designed for making corneal/limbal incision. They are found to be very suitable for making side port.

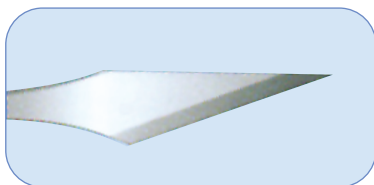

| Product Code | Description | Qty. per box |
|------------------------|--------------|-------------------|
| MIPL/C3/LANCE/SK/15 | 15 Degree | 10 Pcs. (Sterile) |
| MIPL/C3/LANCE/SK/22.5° | 22.50 Degree | 10 Pcs. (Sterile) |
| MIPL/C3/LANCE/SK/30° | 30 Degree | 10 Pcs. (Sterile) |
| MIPL/C3/LANCE/SK/45° | 45 Degree | 10 Pcs. (Sterile) |

MVR Blade

These blades are used for making side port in small incision cataract surgery and for making precise water tight port in vitreoretinal surgery.

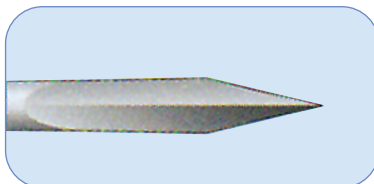

| Product Code | Description | Qty. per box |
|----------------------|-------------|-------------------|
| MIPL/C3/MVR/SK/19G-A | 19 G Angled | 10 pcs. (Sterile) |
| MIPL/C3/MVR/SK/20G-A | 20 G Angled | 10 pcs. (Sterile) |
| MIPL/C3/MVR/SK/23G-A | 23 G Angled | 10 pcs. (Sterile) |
| MIPL/C3/MVR/SK/25G-A | 25 G Angled | 10 pcs. (Sterile) |


Crescent Blade These blades are designed for making scleral tunnel (pocket or groove).



| Product Code | Description | Qty. per box |
|---------------------------|---------------|-------------------|
| MIPL/C3/CRESCENT/SK/2.0-A | 2.0 mm Angled | 10 Pcs. (Sterile) |
| MIPL/C3/CRESCENT/SK/2.5-A | 2.5 mm Angled | 10 Pcs. (Sterile) |

Slit Blade

Precise width keratomes are designed to give incision of required width for insertion of phaco tip, irrigation/aspiration tip.

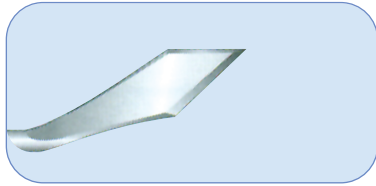

| Product Code | Description | Qty. per box |
|------------------------|-------------|-------------------|
| MIPL/C3/SLIT/SK/1.2A | 1.2 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/1.3-A | 1.3 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/1.4-A | 1.4 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/1.5-A | 1.5 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/1.8-A | 1.8 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/2.0-A | 2.0 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/2.2-A | 2.2 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/2.3-A | 2.3 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/2.4-A | 2.4 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/2.5-A | 2.5 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/2.65-A | 2.65 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/2.75-A | 2.75 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/2.80-A | 2.80 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/3.00-A | 3.00 mm | 10 pcs. (Sterile) |
| MIPL/C3/SLIT/SK/3.20-A | 3.20 mm | 10 pcs. (Sterile) |

Features

- Simple user friendly slide mechanism
- Easy blade visibility during surgery
- Ergonomic design
- Satin finish blades for minimum reflection
- Extremely sharp cutting edge
- Polycarbonate handle
- CE certified

Instructions for Use

The instructions for use and surgical technique are supplied in the form of a leaflet with instructions given in diagram form where appropriate . The international symbols used on the packaging are also explained.

Sterilization

Products are sterilized through Ethylene Oxide (EO) according to a validated cycle.



MADHU INSTRUMENTS PVT LTD

A-260, Okhla Industrial Area, Phase - 1, New Delhi - 110020, India

T: (+91-11) 4270 1030, 31 | E: sales@madhuinstruments.com | www.madhuinstruments.com

Customer Care No. +91 9810866176 (for international queries), +91 7428186111 (for domestic queries)